

Breast Augmentation & Dissectors

Transumbilical



<i>Cat #</i>	<i>Description</i>	<i>Size</i>
BHD	Breast, Hockey stick Dissector	36 cm.
BVD	Breast, "V" Dissector	36 cm.
BDT	Breast, Trocar Dissector	36 cm.
BTDT-S	Breast, Trocar Tube, Small	15mm
BTDT-L	Breast, Trocar Tube, Large	19mm
		Fits BTS & RTS
		Fits BTL & RTL
BTL	Bullet Tip Large	RTL Round Tip Large
RTS	Round Tip Small	BTS Bullet Tip Small

TIRD (Transumbilical Implant Retrieval Device) 5mm diameter, Atraumatic, extra longer grasping device to allow for easy removal of the implant from the umbilicus should the fill tube become separated from the implant. Special "Slide Lock" design handle for perfect grasp each time.

Breast Augmentation & Dissectors

Breast Augmentation



Cat #	Description	Size
BID	Breast, Iconoclast (Straight, Slight curve, Full curve)	12" or 7 1/2"
BCD	Breast, curved trans dissector (Abel Chajchir M.D.)	9"
BPD	Dingman paddle trans dissector, swan-neck	36 cm
DP	Dingman Paddle Slight Curve	36 cm

Breast Augmentation & Dissectors

Agris-Dingman Dissectors



The *Agris-Dingman* stainless steel dissector is a favorite among all breast dissection tools. The *Agris-Dingman* comes standard with a rounded, flattened tip for enhanced dissection.

The tip can be modified with a squared off hook, as used from time to time in specialized procedures.

The *Agris-Dingman* has flattened shafts with incremental markings.

The *Agris-Dingman* is available in pairs only. (Left and Right)

Length: 14" (36 cm.)

Order # ADD

Breast Augmentation & Dissectors

Sub-Glandular Paddle



The K.M.I. Sub-glandular breast paddle is a one piece instrument. The paddle half of the instrument is made of surgical stainless steel and measures 25 cm in length. The edges of the paddle are formed to allow for easy dissection. The handle is made of a high temperature coated aluminum with a contour grip to prevent slipping. The sub-glandular breast paddle is also available in a 2 piece construction for ease of sterilization in short autoclaves.

Order #: SBP

Breast Augmentation & Dissectors

Breast Augmentation Markers

REF: ER-4710 50mm Pattern


REF: ER-4720 45mm Pattern

REF: ER-4630 42mm Pattern

REF: ER-4640 38mm Pattern

Accurately marks amount of tissue for excision. Limbs are semi-malleable for fine adjustments.

Keyhole Markers

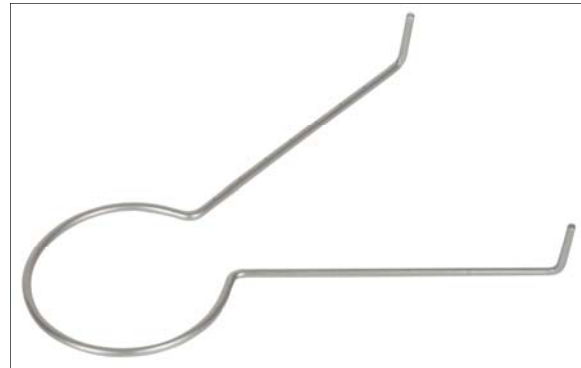

REF: ER-4650 36mm Marker

REF: ER-4660 34mm Marker

REF: ER-4670 50mm Marker

REF: ER-4680 45mm Marker

REF: ER-4690 42mm Marker

REF: ER-4700 38mm Marker

Freeman Areola Marker



Facilitates marking of the areola prior to incision.

Four Ring Marker

REF: EM-4540

Allows for choice of four different diameters; 40mm, 44mm, 48mm and 52mm. Must be used with marking ink.



Breast Augmentation & Dissectors

Gluteus Dissector



The K.M.I. **GLUTEUS** Dissector is designed for intra-muscular dissection. This unique serrated device aides in the dissection of the plane between the Gluteus maximus and the Gluteus medius. The instrument is available in two sizes. The large dissector has a 38 mm diameter disk that is 12.5 mm thick, and the small dissector has a 30 mm diameter disk that is also 12.5 mm thick. The angle of the dissector is ideal for ease of dissection. The solid blue handle is large enough for a comfortable fit in any surgeons hand. The disk and shaft of the dissector are made of 304, full hard, surgical stainless steel. The handle is made of aluminum and is designed for a slip free grip on the device. The standard length of the device is 20 cm ***NOT*** including the handle. Custom lengths are available upon request at no additional cost.

***NOTE:** Device is also available with a screw on/off medical grade delrin handle. This feature allows for sterilization in smaller autoclaves.

Order #:

GDS: Gluteus Dissector Small

GDL: Gluteus Dissector Large